
**Vegetable fats and oils — Determination
of the degradation products of
chlorophylls *a* and *a'* (pheophytins *a*, *a'*
and pyropheophytins)**

*Corps gras d'origine végétale — Détermination des produits de
décomposition des chlorophylles *a* et *a'* (phéophytines *a*, *a'* et
pyrophéophytines)*



PDF disclaimer

This PDF file may contain embedded typefaces. In accordance with Adobe's licensing policy, this file may be printed or viewed but shall not be edited unless the typefaces which are embedded are licensed to and installed on the computer performing the editing. In downloading this file, parties accept therein the responsibility of not infringing Adobe's licensing policy. The ISO Central Secretariat accepts no liability in this area.

Adobe is a trademark of Adobe Systems Incorporated.

Details of the software products used to create this PDF file can be found in the General Info relative to the file; the PDF-creation parameters were optimized for printing. Every care has been taken to ensure that the file is suitable for use by ISO member bodies. In the unlikely event that a problem relating to it is found, please inform the Central Secretariat at the address given below.



COPYRIGHT PROTECTED DOCUMENT

© ISO 2009

All rights reserved. Unless otherwise specified, no part of this publication may be reproduced or utilized in any form or by any means, electronic or mechanical, including photocopying and microfilm, without permission in writing from either ISO at the address below or ISO's member body in the country of the requester.

ISO copyright office
Case postale 56 • CH-1211 Geneva 20
Tel. + 41 22 749 01 11
Fax + 41 22 749 09 47
E-mail copyright@iso.org
Web www.iso.org

Published in Switzerland

Foreword

ISO (the International Organization for Standardization) is a worldwide federation of national standards bodies (ISO member bodies). The work of preparing International Standards is normally carried out through ISO technical committees. Each member body interested in a subject for which a technical committee has been established has the right to be represented on that committee. International organizations, governmental and non-governmental, in liaison with ISO, also take part in the work. ISO collaborates closely with the International Electrotechnical Commission (IEC) on all matters of electrotechnical standardization.

International Standards are drafted in accordance with the rules given in the ISO/IEC Directives, Part 2.

The main task of technical committees is to prepare International Standards. Draft International Standards adopted by the technical committees are circulated to the member bodies for voting. Publication as an International Standard requires approval by at least 75 % of the member bodies casting a vote.

Attention is drawn to the possibility that some of the elements of this document may be the subject of patent rights. ISO shall not be held responsible for identifying any or all such patent rights.

ISO 29841 was prepared by Technical Committee ISO/TC 34, *Food products*, Subcommittee SC 11, *Animal and vegetable fats and oils*.

Vegetable fats and oils — Determination of the degradation products of chlorophylls *a* and *a'* (pheophytins *a*, *a'* and pyropheophytins)

1 Scope

This International Standard specifies a procedure for the determination of the degradation products pheophytin *a*, *a'* and pyropheophytin *a* of chlorophylls. The method is applicable to vegetable fats and oils only.

2 Normative references

The following referenced documents are indispensable for the application of this document. For dated references, only the edition cited applies. For undated references, the latest edition of the referenced document (including any amendments) applies.

ISO 661, *Animal and vegetable fats and oils — Preparation of test sample*

3 Terms and definitions

For the purposes of this document, the following terms and definitions apply.

3.1

pyropheophytin *a* (content in chlorophyll degradation products)

mass of pyropheophytin *a* divided by the sum of the masses of pheophytin *a*, *a'* and pyropheophytin *a*

4 Principle

Separation of the pigments (chlorophylls, pheophytins *a* and *a'*, pyropheophytin *a*) using a miniaturized column chromatography on a silica gel column from the major part of the lipids (triglycerides). The eluate is analysed by RP18-HPLC, and the separated components are monitored at 410 nm using a photometric detector.

The analysis of some seed oils (e.g. rapeseed oil) requires a higher sensitivity which can be achieved using fluorescence detection with an excitation wavelength of 430 nm and an emission wavelength of 670 nm.

A direct HPLC determination with fluorescence detection is also possible and described in Annex C.

5 Reagents

WARNING — Attention is drawn to the regulations which specify the handling of hazardous substances. Technical, organizational and personal safety measures shall be followed.

During the analysis, unless otherwise stated, use only reagents of recognized analytical grade.

5.1 Chlorophyll a.

NOTE Due to the lack of pheophytin and pyropheophytin reference materials, chlorophyll a is used to produce these degradation products. Pheophytins are obtained from a chlorophyll a solution by acidification; pyropheophytins by additional heat ($> 100\text{ }^{\circ}\text{C}$) (Reference [4]).

5.2 Acetone.

5.3 Methanol.

5.4 Diethyl ether.

5.5 Petroleum ether, boiling point range $40\text{ }^{\circ}\text{C}$ to $60\text{ }^{\circ}\text{C}$.

5.6 *n*-Heptane or, if not available, *n*-hexane.

CAUTION — Owing to its toxicity, avoid using *n*-hexane if *n*-heptane is available.

5.7 Petroleum ether and diethyl ether solvent mixture, volume fractions: $\varphi_{\text{PetE}} = 90\text{ ml}/100\text{ ml}$; $\varphi_{\text{Et}_2\text{O}} = 10\text{ ml}/100\text{ ml}$.

6 Apparatus

Usual laboratory equipment, and in particular the following.

6.1 HPLC setup.

6.1.1 HPLC system, consisting of a pump, a sample-injecting device ($20\text{ }\mu\text{l}$ loop) or auto sampler, either a photometric detector for measurements at 410 nm or a fluorescence detector (excitation wavelength, $\lambda_{\text{ex}} = 430\text{ nm}$ and emission wavelength, $\lambda_{\text{em}} = 670\text{ nm}$), and an integration system.

6.1.2 HPLC column, 250 mm length, $4,0\text{ mm}$ or $4,6\text{ mm}$ internal diameter, filled with reverse-phase type RP18 particles, size $5\text{ }\mu\text{m}$.

The following conditions have been found to be suitable:

Stationary phase: Partisil ODS 3, $5\text{ }\mu\text{m}$, $250\text{ mm} \times 4,6\text{ mm}$ ¹⁾;

Mobile phase: Water, methanol, and acetone (6.1.3);

Flow rate: $1,0\text{ ml}/\text{min}$.

6.1.3 Mobile phase, water, methanol, and acetone mixture, volume fractions: $\varphi_{\text{H}_2\text{O}} = 4\text{ ml}/100\text{ ml}$; $\varphi_{\text{MeOH}} = 36\text{ ml}/100\text{ ml}$; $\varphi_{\text{Me}_2\text{C=O}} = 60\text{ ml}/100\text{ ml}$. All solvents shall be HPLC grade.

6.2 Taper-shaped flask, of capacity 10 ml or 25 ml .

6.3 Beakers, of different capacities.

6.4 Rotary evaporator, with water bath.

6.5 Auto sampler vials, of suitable capacity.

1) Example, available commercially, of a C18 reverse phase with a $10,5\%$ mass fraction carbon load, end-capped for deactivation of silanols to minimize the need for ion suppression or ion pairing agents. This information is given for the convenience of users of this International Standard and does not constitute an endorsement by ISO of this product.

6.6 Silica cartridge, 1 000 mg/6 ml, 55 µm, 700 nm or **diol cartridge**, 3 ml ²⁾.

NOTE Home-made silica mini-columns can also be used for the separation.

7 Sample

7.1 Sampling

Sampling is not part of the method specified in this International Standard. A recommended sampling method is given in ISO 5555 [1].

It is important that the laboratory receive a sample which is truly representative and has not been damaged or changed during transport or storage.

7.2 Preparation of the test sample

Prepare the test sample in accordance with ISO 661.

IMPORTANT — Test samples shall be stored cold and protected from light in glass vessels. After clean-up, the test solutions shall be analysed *immediately* to prevent degradation and oxidation. Avoid contact of test samples with plastics.

8 Procedure

8.1 Weigh about 300 mg of the test sample (7.2) into a small beaker (6.3), dissolve it in 1 ml of *n*-heptane or *n*-hexane (5.6) and pour the solution on to the silica cartridge (silica gel column) (6.6). Rinse the beaker twice with 1 ml portions of petroleum ether (5.5), pouring these washings on to the cartridge (column) also.

8.2 As soon as the solvent has drained to the top of the cartridge (column), elute the non-polar substances twice with 5 ml petroleum ether/diethyl ether solvent mixture (5.7) each time.

8.3 Elute the pheophytin fraction twice with 5 ml acetone (5.2) and collect this fraction in a taper-shaped flask (6.2) protected from light.

8.4 Evaporate the solvent to dryness on a rotary evaporator (6.4) at max 20 °C. Dissolve the residue in 200 µl acetone. Introduce this solution immediately into the HPLC setup (6.1).

IMPORTANT — Pheophytins are very unstable in light and oxidize easily.

8.5 Typical chromatograms run under the column conditions listed in 6.1.2 are shown in Figures A.1 and A.2.

9 Results of the determination

9.1 Peak identification

The chromatographic pattern of the determination may show the peaks of the pigments chlorophyll *a*, chlorophyll *b*, the corresponding pheophytins and pyropheophytins. Reference should be made to the typical chromatograms presented in Figures A.1 and A.2 as the other pigments are not available as reference

2) Strata Si-1 (silica cartridge) and Supelclean LC-Diol SPE Tubes – Supelco 57016 (diol cartridge) are examples of suitable products available commercially. This information is given for the convenience of users of this International Standard and does not constitute an endorsement by ISO of these products.

materials due to their instability. To identify pheophytin *a*, pheophytin *a'* and pyropheophytin *a*, see Figures A.1 and A.2.

9.2 Calculation of the pyropheophytin *a* content

Use the peak areas to calculate the relative proportions of the analytes in the sample solution. For the purpose of this method, it is assumed that the response factors of all pigments are equal. The pyropheophytin *a* content as a percentage mass fraction, w_{PPPa} , is calculated on the basis of a peak area ratio as follows:

$$w_{\text{PPPa}} = \frac{A_{\text{PPPa}} \times 100}{(A_{\text{PPPa}} + \sum A_{\text{PP}})}$$

where

A_{PPPa} is the peak area of pyropheophytin *a*;

$\sum A_{\text{PP}}$ is the sum of the peak areas of pheophytin *a* and *a'*.

Report the result to one decimal place.

10 Precision of the method

10.1 Interlaboratory test

Details of an interlaboratory test on the precision of the method are summarized in Annex B. The values derived from this interlaboratory test may not be applicable to concentration ranges and matrices other than those given.

10.2 Repeatability

The absolute difference between two independent single test results, obtained with the same method on identical test material in the same laboratory by the same operator using the same equipment within a short interval of time, will in not more than 5 % of cases exceed the value of *r* given in Table B.1.

10.3 Reproducibility

The absolute difference between two single test results, obtained with the same method on identical test material in different laboratories by different operators using different equipment, will in not more than 5 % of cases exceed the value of *R* given in Table B.1.

11 Test report


The test report shall contain at least the following information:

- a) all information necessary for the complete identification of the sample;
- b) the sampling method used, if known;
- c) the test method used, with reference to this International Standard;

- d) all operating details not specified in this International Standard, or regarded as optional, together with details of any incidents that may have occurred when performing the method, which may have influenced the test result(s);
- e) the test result obtained;
- f) if the repeatability has been checked, the final quoted result obtained.

Annex A (informative)


Example of a typical chromatogram

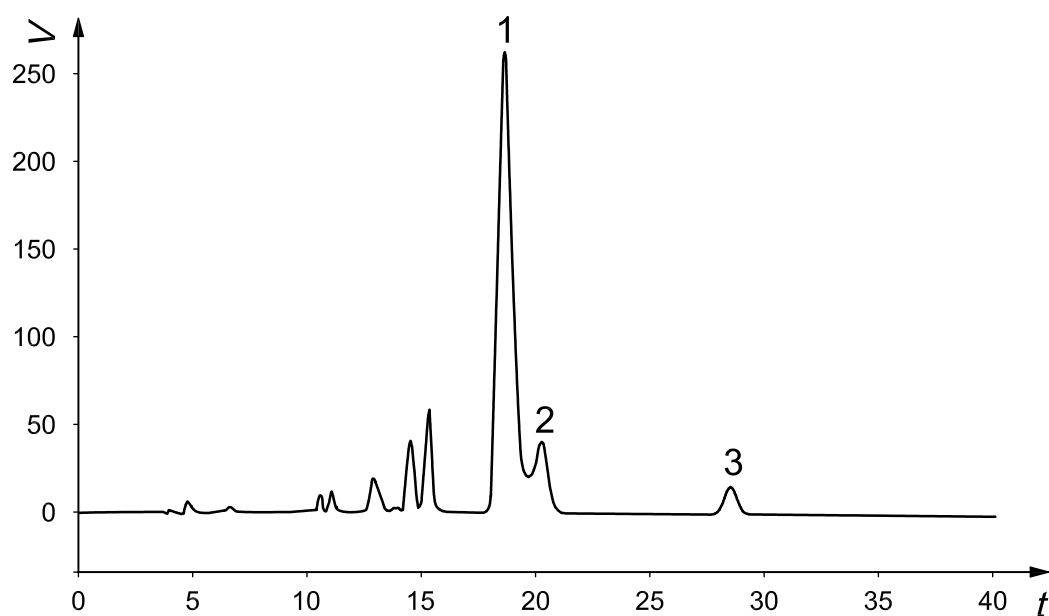

Key

- 1 pheophytin *a*
- 2 pheophytin *a'*
- 3 pyropheophytin
- t* elution time, min
- V* potential difference, volts

Figure A.1 — HPLC chromatogram of an extra virgin olive oil with 5 % mass fraction pyropheophytin content in chlorophyll degradation products

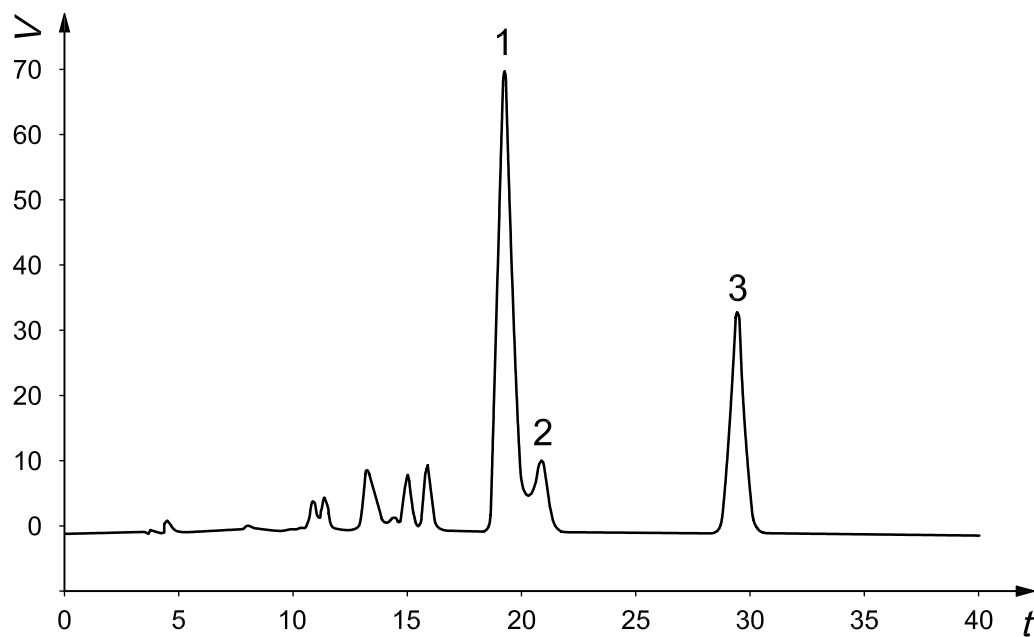**Key**1 pheophytin *a*2 pheophytin *a'*

3 pyropheophytin

t elution time, min*V* potential difference, volts

Figure A.2 — HPLC chromatogram of an extra virgin olive oil with 30 % mass fraction pyropheophytin content in chlorophyll degradation products

Annex B (informative)

Results of an interlaboratory test

The precision of the method is the result of an interlaboratory study organized by the Joint Committee for the Analysis of Fats, Oils, Fat Products, Related Products and Raw Materials (DIN/DGF) in which 16 laboratories from France, Germany, and Italy participated. The study was carried out in 2005 on five samples, and precision data were calculated in accordance with ISO 5725-1 [2] and ISO 5725-2 [3].

Table B.1 — Summary of statistical results of olive oils

Sample	A	B	C	D	E
Number of participating laboratories, N	16	16	16	16	16
Number of laboratories retained after eliminating outliers, n	12	12	14	15	15
Number of individual test results of all laboratories on each sample, $\sum_{i=1}^n n_{w,i}$	24	24	28	30	30
Mean value, \bar{w}_{PPP2} , %	28,86	34,01	5,63	6,26	84,75
Repeatability standard deviation, s_r , %	0,53	0,61	0,59	0,36	1,09
Repeatability coefficient of variation, $CV(r)$, %	1,8	1,8	10,4	5,8	1,3
Repeatability limit, r ($s_r \times 2,8$), %	1,48	1,71	1,65	1,01	3,05
Reproducibility standard deviation, s_R , %	1,78	2,08	0,79	1,29	3,06
Reproducibility coefficient of variation, $CV(R)$, %	6,2	6,1	14,1	20,7	3,6
Reproducibility limit, R ($s_R \times 2,8$), %	4,98	5,82	2,21	3,61	8,56

Annex C (informative)

Direct HPLC method with fluorescence detection

C.1 Principle

This annex describes an alternative method. The oils are submitted to HPLC without any prior treatment.

C.2 Solvents for HPLC

C.2.1 Water, for HPLC.

C.2.2 Methanol, for HPLC.

C.2.3 Acetone, for HPLC.

C.3 Apparatus

C.3.1 Binary HPLC system, consisting of a pump, a sample-injecting device (20 µl loop) or auto sampler, a fluorescence detector ($\lambda_{\text{ex}} = 430 \text{ nm}$ and $\lambda_{\text{em}} = 670 \text{ nm}$), and an integration system.

C.3.2 HPLC column, 250 mm length, 4,0 mm or 4,6 mm internal diameter, filled with reverse-phase type RP18 particles, size 5 µm.


C.4 Procedure

C.4.1 Transfer the oils into an amber HPLC vial.

C.4.2 Introduce the test portion into the HPLC apparatus (Clause C.3). The conditions listed in Table C.1 have been found to be suitable and to give a typical chromatogram, such as that shown in Figure C.1.

Table C.1 — Elution programme with 1,0 ml/min

Time (min)	Solvent A	Solvent B
	Volume fractions: $\varphi_{\text{H}_2\text{O}} = 4 \%$; $\varphi_{\text{MeOH}} = 76 \%$; $\varphi_{\text{Me}_2\text{C}=\text{O}} = 20 \%$	Volume fractions: $\varphi_{\text{H}_2\text{O}} = 4 \%$; $\varphi_{\text{MeOH}} = 36 \%$; $\varphi_{\text{Me}_2\text{C}=\text{O}} = 60 \%$
0	100 %	0 %
7	0 %	100 %
50	0 %	100 %
50,1	100 %	0 %
60	100 %	0 %

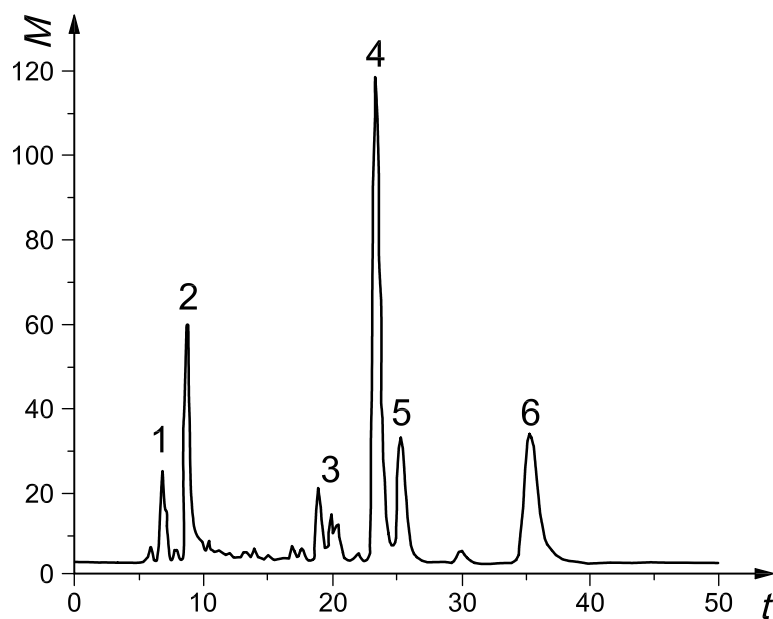

Key

- 1 pheophorbide
- 2 methylpheophorbide
- 3 pheophytin *b*
- 4 pheophytin *a*
- 5 pheophytin *a'*
- 6 pyropheophytin

M fluorescence detection, lux

t elution time, min

Figure C.1 — HPLC chromatogram of a rapeseed oil with pyropheophytin

Bibliography

- [1] ISO 5555, *Animal and vegetable fats and oils — Sampling*
- [2] ISO 5725-1, *Accuracy (trueness and precision) of measurement methods and results — Part 1: General principles and definitions*
- [3] ISO 5725-2, *Accuracy (trueness and precision) of measurement methods and results — Part 2: Basic method for the determination of repeatability and reproducibility of a standard measurement method*
- [4] AITZETMÜLLER, K. Chlorophyll-Abbauprodukte in pflanzlichen Ölen [Chlorophyll degradation products in vegetable oils]. *Fat Sci. Technol.* 1989, **91**, pp. 99-105; AITZETMÜLLER, K. Analysis of chlorophyll degradation products in fats and oilseeds. In: *Actes du congrès international «Chevreul» pour l'étude des corps gras* [Proceedings of the Chevreul International Congress], Angers, 1989-06-06/10, p. 177-84. Association française pour l'Étude des Corps gras, Neuilly sur Seine

