

---

---

**Cheese — Determination of nisin A  
content by LC-MS and LC-MS-MS**

*Fromage — Détermination de la teneur en nisine A par CL-SM et  
CL-SM-SM*



Reference numbers  
ISO/TS 27106:2009(E)  
IDF/RM 217:2009(E)

© ISO and IDF 2009

**PDF disclaimer**

This PDF file may contain embedded typefaces. In accordance with Adobe's licensing policy, this file may be printed or viewed but shall not be edited unless the typefaces which are embedded are licensed to and installed on the computer performing the editing. In downloading this file, parties accept therein the responsibility of not infringing Adobe's licensing policy. Neither the ISO Central Secretariat nor the IDF accepts any liability in this area.

Adobe is a trademark of Adobe Systems Incorporated.

Details of the software products used to create this PDF file can be found in the General Info relative to the file; the PDF-creation parameters were optimized for printing. Every care has been taken to ensure that the file is suitable for use by ISO member bodies and IDF national committees. In the unlikely event that a problem relating to it is found, please inform the ISO Central Secretariat at the address given below.



**COPYRIGHT PROTECTED DOCUMENT**

© ISO and IDF 2009

All rights reserved. Unless otherwise specified, no part of this publication may be reproduced or utilized in any form or by any means, electronic or mechanical, including photocopying and microfilm, without permission in writing from either ISO or IDF at the respective address below.

ISO copyright office  
Case postale 56 • CH-1211 Geneva 20  
Tel. + 41 22 749 01 11  
Fax + 41 22 749 09 47  
E-mail [copyright@iso.org](mailto:copyright@iso.org)  
Web [www.iso.org](http://www.iso.org)

International Dairy Federation  
Diamant Building • Boulevard Auguste Reyers 80 • B-1030 Brussels  
Tel. + 32 2 733 98 88  
Fax + 32 2 733 04 13  
E-mail [info@fil-idf.org](mailto:info@fil-idf.org)  
Web [www.fil-idf.org](http://www.fil-idf.org)

Published in Switzerland

## Foreword

**ISO (the International Organization for Standardization)** is a worldwide federation of national standards bodies (ISO member bodies). The work of preparing International Standards is normally carried out through ISO technical committees. Each member body interested in a subject for which a technical committee has been established has the right to be represented on that committee. International organizations, governmental and non-governmental, in liaison with ISO, also take part in the work. ISO collaborates closely with the International Electrotechnical Commission (IEC) on all matters of electrotechnical standardization.

International Standards are drafted in accordance with the rules given in the ISO/IEC Directives, Part 2.

The main task of technical committees is to prepare International Standards. Draft International Standards adopted by the technical committees are circulated to the member bodies for voting. Publication as an International Standard requires approval by at least 75 % of the member bodies casting a vote.

In other circumstances, particularly when there is an urgent market requirement for such documents, a technical committee may decide to publish other types of document:

- an ISO Publicly Available Specification (ISO/PAS) represents an agreement between technical experts in an ISO working group and is accepted for publication if it is approved by more than 50 % of the members of the parent committee casting a vote;
- an ISO Technical Specification (ISO/TS) represents an agreement between the members of a technical committee and is accepted for publication if it is approved by 2/3 of the members of the committee casting a vote.

An ISO/PAS or ISO/TS is reviewed after three years in order to decide whether it will be confirmed for a further three years, revised to become an International Standard, or withdrawn. If the ISO/PAS or ISO/TS is confirmed, it is reviewed again after a further three years, at which time it must either be transformed into an International Standard or be withdrawn.

Attention is drawn to the possibility that some of the elements of this document may be the subject of patent rights. ISO shall not be held responsible for identifying any or all such patent rights.

ISO/TS 27106|IDF/RM 217 was prepared by Technical Committee ISO/TC 34, *Food products*, Subcommittee SC 5, *Milk and milk products*, and the International Dairy Federation (IDF). It is being published jointly by ISO and IDF.

## Foreword

**IDF (the International Dairy Federation)** is a non-profit organization representing the dairy sector worldwide. IDF membership comprises National Committees in every member country as well as regional dairy associations having signed a formal agreement on cooperation with IDF. All members of IDF have the right to be represented on the IDF Standing Committees carrying out the technical work. IDF collaborates with ISO in the development of standard methods of analysis and sampling for milk and milk products.

The main task of Standing Committees is to prepare International Standards. Draft International Standards adopted by the Standing Committees are circulated to the National Committees for endorsement prior to publication as an International Standard. Publication as an International Standard requires approval by at least 50% of IDF National Committees casting a vote.

In other circumstances, particularly when there is an urgent market requirement for such documents, a Standing Committee may decide to publish an other type of normative document which is called by IDF: *Reviewed method*. Such a method represents an agreement between the members of a Standing Committee and is accepted for publication if it is approved by at least 50 % of the committee members casting a vote. A *Reviewed method* is equal to an ISO/PAS or ISO/TS and will, therefore, also be published jointly under ISO conditions.

Attention is drawn to the possibility that some of the elements of this document may be the subject of patent rights. IDF shall not be held responsible for identifying any or all such patent rights.

ISO/TS 27106|IDF/RM 217 was prepared by the International Dairy Federation (IDF) and Technical Committee ISO/TC 34, *Food products*, Subcommittee SC 5, *Milk and milk products*. It is being published jointly by IDF and ISO.

All work was carried out by the Joint ISO-IDF Action Team on *Food additives and vitamins* of the Standing Committee on *Analytical methods for additives and contaminants* under the aegis of its project leader, Mr. T. Berger (CH).

# Cheese — Determination of nisin A content by LC-MS and LC-MS-MS

## 1 Scope

This Technical Specification specifies a method for the quantitative determination of the nisin A content in cheese.

The method is suitable for measuring low levels of nisin A with a quantification limit of 1 mg/kg.

**NOTE** Nisin is a peptide produced by some bacteria (e.g. *Lactococcus lactis* subsp. *Lactis*) inhibiting or destroying other microorganisms. It is widely used as a natural preservative for foods, e.g. vegetables, cheese, meat, and cacao. In cheese making, nisin is used to prevent late blowing. Its use is restricted to maximum levels in the final product. Nisin appears in two forms, nisin A and nisin Z, which differ in one amino acid. This method is applicable to the determination of nisin A only.

## 2 Terms and definitions

For the purposes of this document, the following terms and definitions apply.

### 2.1

#### **nisin A content**

mass fraction of substance determined by the procedure specified in this Technical Specification

**NOTE** The nisin A content is expressed in milligrams per kilogram.

## 3 Principle

The sample is grated and extracted with dilute formic acid at 80 °C. After ultracentrifugation, interfering proteins are separated by means of filtration through an ultrafiltration (UF) membrane. In the purified extract, nisin A is separated using a polymeric stationary phase and detected using mass spectrometry (with MS-MS as an option).

## 4 Reagents and reference substances

Use only reagents of recognized analytical grade and distilled water or water of at least equivalent purity, unless otherwise specified.

**4.1 Bovine serum albumin (BSA) stock solution.** Dissolve 10 mg of BSA (purity > 96 % mass fraction), in 10 ml water.

**4.2 Bovine serum albumin (BSA) buffer solution.** Mix 80 ml water with 20 ml of acetonitrile (4.6), 0,5 ml of formic acid (4.3), 0,01 ml of trifluoroacetic acid (4.5) and 1 ml of BSA stock solution (4.1).

**4.3 Formic acid** (HCOOH).

**4.4 Formic acid solution,**  $\rho_{\text{HCOOH}} = 5 \text{ g/l}$ . Pipette 0,41 ml of formic acid (4.3) into a 100 ml one-mark volumetric flask (5.12). Make up to the mark with water and mix.

- 4.5 **Trifluoroacetic acid** ( $\text{CF}_3\text{COOH}$ ).
- 4.6 **Acetonitrile** ( $\text{CH}_3\text{CN}$ ), “pure”.
- 4.7 **Methanol** ( $\text{CH}_3\text{OH}$ ).
- 4.8 **Nisin A**, of purity > 95 % mass fraction<sup>1)</sup>.

## 5 Apparatus

Usual laboratory equipment and, in particular, the following.

- 5.1 **Laboratory centrifuge**, capable of producing a radial acceleration of at least 3 000g.
- 5.2 **Ultracentrifuge**, capable of producing a radial acceleration of 20 800g.
- 5.3 **Ultrafiltration membrane**, of pore size 30 kD<sup>2)</sup>.
- 5.4 **Membrane filter**, of pore size 0,22  $\mu\text{m}$ <sup>3)</sup>.
- 5.5 **Balance**, capable of weighing to the nearest 10 mg, with a readability of 1 mg.
- 5.6 **Analytical balance**, capable of weighing to the nearest 0,1 mg, with a readability of 0,01 mg.
- 5.7 **Water bath**, capable of shaking and maintaining a temperature of  $80\text{ }^\circ\text{C} \pm 2\text{ }^\circ\text{C}$ .
- 5.8 **Ultrasonic bath**, capable of shaking and maintaining a temperature of  $80\text{ }^\circ\text{C} \pm 2\text{ }^\circ\text{C}$ .
- 5.9 **Cheese grater**, with openings of size approximately 2 mm.
- 5.10 **LC-MS equipment**.
  - 5.10.1 **Elution gradient pumping system**, capable of operating at 0,25 ml/min.
  - 5.10.2 **Manual or automatic injector**, capable of injecting volumes of 5  $\mu\text{l}$ .
  - 5.10.3 **Column heater**, capable of maintaining a column temperature of  $40\text{ }^\circ\text{C} \pm 2\text{ }^\circ\text{C}$ .
  - 5.10.4 **Column, reversed phase**, PLRP-S, 300 Å<sup>4)</sup>, 3  $\mu\text{m}$ , 150 mm  $\times$  2 mm.
  - 5.10.5 **Mass spectrometer detector**, capable of operating in ion mode ESI+ at  $m/z$  839,6.
- 5.11 **LC-MS-MS equipment** (optional).

---

1) Ambicin N is an example of a suitable product available commercially. This information is given for the convenience of users of this document and does not constitute an endorsement by ISO or IDF of this product.

2) Millipore Centricon YM-30 is an example of a suitable product available commercially. This information is given for the convenience of users of this document and does not constitute an endorsement by ISO or IDF of this product.

3) Millipore Millex-GV PVDF 0,22  $\mu\text{m}$  is an example of a suitable product available commercially. This information is given for the convenience of users of this document and does not constitute an endorsement by ISO or IDF of this product.

4) PLRP-S, 300 Å is the trade name of a product supplied by Polymer Lab Ltd. This information is given for the convenience of users of this document and does not constitute an endorsement by ISO or IDF of the product named. Equivalent products may be used if they can be shown to lead to the same results.

**5.11.1 Elution gradient pumping system**, capable of delivering a flow rate of 0,2 ml/min.

**5.11.2 Manual or automatic injector**, capable of injecting volumes of 10 µl.

**5.11.3 Column heater**, capable of maintaining a column temperature of 30 °C ± 2 °C.

**5.11.4 Column for reversed phase chromatography**, PLRP-S, 300 Å<sup>4</sup>, 3 µm, 150 × 2 mm.

**5.11.5 Mass spectrometer detector**, capable of operating in ion mode ESI+ MS-MS at *m/z* 672/672, 672/811, 672/649, 840/840.

**5.12 One-mark volumetric flasks**, of capacity 100 ml, ISO 1042<sup>[2]</sup> class A.

## 6 Sampling

A representative sample should have been sent to the laboratory. It should not have been damaged or changed during transport or storage.

Sampling is not part of the method specified in this Technical Specification. A recommended sampling method is given in ISO 707 | IDF 50<sup>[1]</sup>.

## 7 Procedure

### 7.1 Preparation of the nisin A standard solution

**7.1.1 Nisin A standard stock solution**,  $\rho_{nA} = 100$  mg/l.

Weigh, to the nearest 0,01 mg, 10,00 mg nisin A (4.8) into a 100 ml one-mark volumetric flask (5.12). Make up to the mark with formic acid solution (4.4) and mix.

Prepare the nisin A standard stock solution daily.

**7.1.2 Nisin A standard working solution**,  $\rho_{nA} = 300$  µg/l.

Pipette 300 µl of nisin A standard stock solution (7.1.1) into a 100 ml one-mark volumetric flask (5.12). Make up to the mark with BSA buffer solution (4.2) and mix. The standard working solution thus obtained contains 300 µg of nisin A per litre.

Prepare the nisin A standard working solution daily.

### 7.2 Extraction of the test portion

Before weighing, grate the test sample with the cheese grater (5.9).

Weigh, to the nearest 0,01 g, 2,50 g of test sample into a 100 ml one-mark volumetric flask (5.12). Add 70 ml of water and 0,5 ml of formic acid (4.3) and mix.

Put the flask either in the water bath (5.7) at 80 °C and shake it for 30 min or in the ultrasonic bath (5.8) at 80 °C and shake it for 10 min. After cooling to room temperature, make up to the mark with water and mix.

NOTE Soft cheese can be grated after freezing.

### 7.3 Filtration of the test portion with UF membrane

Pipette approximately 1,5 ml of the extract into a 1,5 ml tube (e.g. Eppendorf) and centrifuge using the ultracentrifuge (5.2) at 20 800g for 10 min.

Determine the tare of the collecting container for the ultrafiltration membrane (5.3) on the analytical balance (5.6). Place the collecting container on the ultrafiltration membrane and set the analytical balance to zero.

Filter 0,6 ml to 0,7 ml of the centrifuged extract through a membrane filter (5.4) into the tared ultra filtration membrane.

Weigh the amount of the extract on the ultrafiltration membrane using the analytical balance. Centrifuge the ultrafiltration membrane in the laboratory centrifuge (5.1) at 3 000g for 45 min. Determine the gross mass of the filter container.

Supplement the net mass of filtrate with water to the original amount weighed with a correction for the peptide distribution in the ultrafiltration membrane.

NOTE Generally, the amount of water added is approximately 15 % of the extract used for centrifugation. Consult the information provided by the membrane manufacturer.

Transfer the filtrate thus obtained into an HPLC vial and measure.

### 7.4 LC-MS and LC-MS-MS determination

#### 7.4.1 Elution solvents for LC-MS

Use the following elution solvents:

Eluant A: pipette 2,5 ml of formic acid (4.3) and 0,05 ml of trifluoroacetic acid (4.5) into 500 ml water.

Eluant B: mix 350 ml of acetonitrile (4.6) with 150 ml water, 2,5 ml of formic acid (4.3) and 0,05 ml of trifluoroacetic acid (4.5).

#### 7.4.2 Elution solvents for LC-MS-MS

Optionally, the following elution solvents can be used if LC-MS-MS is chosen:

Eluant C: mix 50 ml acetonitrile (4.6) with 450 ml water and 2,5 ml formic acid (4.3).

Eluant D: mix 400 ml acetonitrile (4.6) with 100 ml water and 2,5 ml formic acid (4.3).




#### 7.4.3 LC-MS and LC-MS-MS conditions

Table 1 — Preferred conditions

Conditions	LC-MS			LC-MS-MS		
Injection volume, µl	5			10		
Column	PLRP-S <sup>4)</sup> 150 mm × 2 mm, 300 Å, 3 µm					
Column temperature, °C	40			30		
Mass detector	Ion mode ESI+ Sample 500 °C Needle 3,5 kV			Ion mode ESI+ Nitrogen 8 l/min 350 °C		
Detection	<i>m/z</i> 839,6 <sup>a</sup> Span 0,5 Dwell time 1,0 min Cone 60 V			<i>m/z</i> Q1/ <i>m/z</i> Q3 672/672 <sup>b</sup> 672/811 <sup>c</sup> (most selective transition) 672/649 840/840 <sup>b</sup>		
Flow rate, ml/min	0,25			0,20		
Nisin A retention time, min	~13			~5,7		
Gradient <sup>d</sup>	Min.	%A	%B	Min.	%C	%D
	0	70	30	0	85	15
	13	50	50	1	85	15
	13,1	0	100	4	65	35
	17	0	100	6	65	35
	17,1	70	30	9	10	90
	23	70	30	12	10	90
				12,1	85	15
				16	85	15
On the whole, the MS or MS-MS parameters are critical and shall be carefully controlled to achieve sensitivity and repeatability because of different analytical factors. Optimize the instrument to produce the highest signal for nisin.						
<sup>a</sup> Higher sensitivities for ions <i>m/z</i> 1 118,9+3 ([M+3H] <sup>3+</sup> ) or <i>m/z</i> 1 677,7+2 ([M+2H] <sup>2+</sup> ) were observed in other laboratories and could be used alternatively.						
<sup>b</sup> Identical <i>m/z</i> signal, equipment for MS-MS has to be optimized.						
<sup>c</sup> Due to a lower ion charge, a higher <i>m/z</i> signal occurs in Q3.						
<sup>d</sup> The elution gradient might require slight modification in order to achieve the resolution shown in Figure 1.						

#### 7.4.4 LC-MS calibration

An example of a chromatogram of the nisin A standard solution is given in Figure 1. The chromatogram was obtained by measuring the nisin A standard working solution (7.1.2) three times. The example shows that calibration is linear for the whole measuring range (see Figure 2). The equipment linearity, however, should be checked regularly.

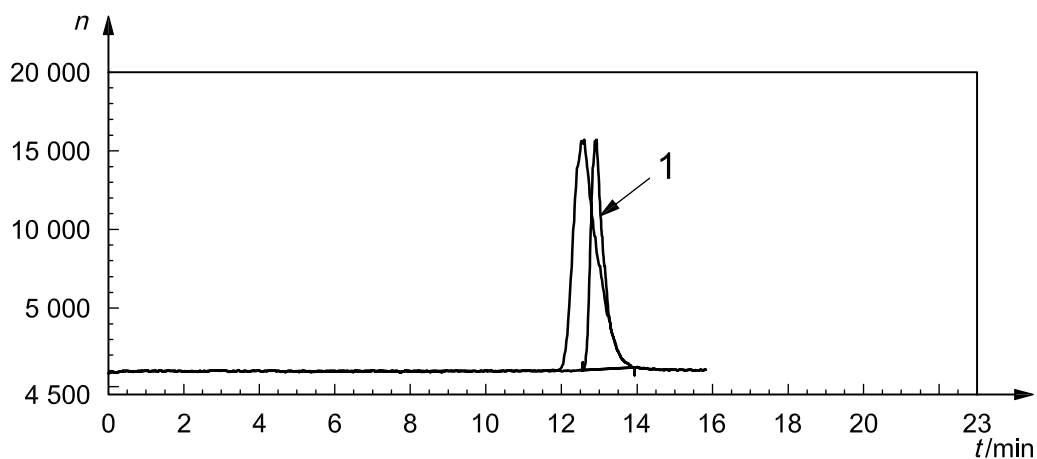

**Key**


$n$  counts

$t$  time

1 nisin A

SIM\_01 60,0 positive; 839,34 to 839,84  $m/z$

**Figure 1 — Nisin A standard working solution with LC-MS**

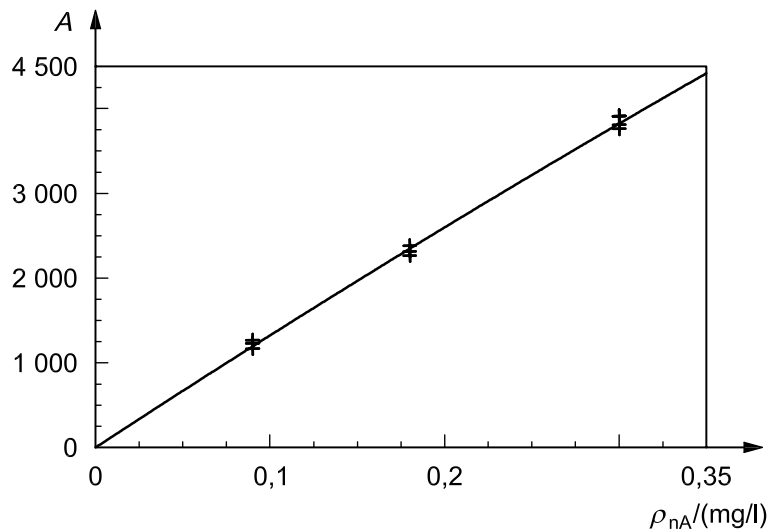

**Key**

$A$  area (counts times time, in minutes)


$\rho_{nA}$  nisin A mass concentration

[nisin + 4H]<sup>4+</sup> external SIM\_01

**Figure 2 — Nisin A calibration at  $m/z$  839,6 with LC-MS**

#### 7.4.5 LC-MS-MS calibration

The transition  $m/z$  672 to 811 is most selective for nisin A and, therefore, was used for the calibration. Figure 3 shows the result of linear regression obtained by using four concentrations.

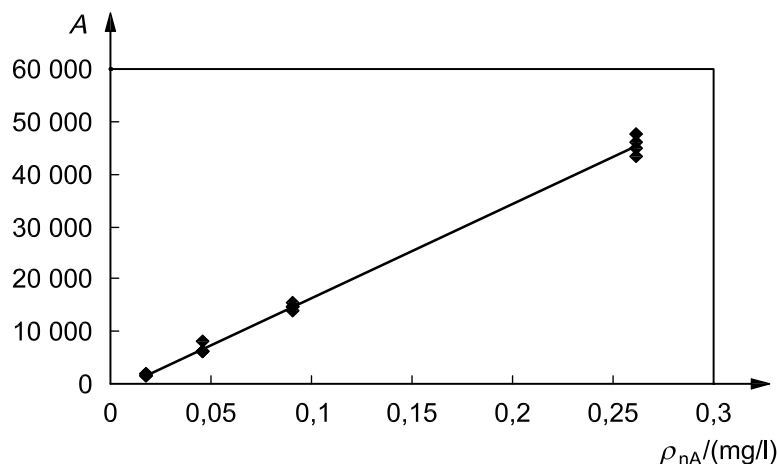

#### Key

$A$  area counts

$\rho_{nA}$  nisin A mass concentration

$$A = 180\,342\rho_{nA} - 1\,621,2$$

$$r^2_{\rho_{nA}A} = 0,996\,7$$

Figure 3 — Nisin A calibration with LC-MS-MS

#### 7.4.6 Determination of the nisin A content in the sample

Always perform a single-point calibration procedure in every series of analyses. If necessary, recalibrate the system and check the reagent blank.

### 8 Calculation and expression of results

#### 8.1 Calculation

Calculate the nisin A content,  $w_{nA}$ , of the test sample, expressed in milligrams of nisin A per kilogram of test sample, by using the following equation:

$$w_{nA} = \frac{H_t \rho_{nA, s} V_t}{H_s m_t}$$

where

$\rho_{nA, s}$  is the concentration, in milligrams per litre, of the standard working solution;

$H_t$  is the numerical value of the peak height or area of the test solution;

$H_s$  is the numerical value of the peak height or area of the standard working solution;

$m_t$  is the mass, in grams, of cheese used for the preparation of the test solution;

$V_t$  is the volume, in millilitres, of the test solution.

## 8.2 Expression of results

Express the results to one decimal place. Express test results of below 1 mg/kg as “less than 1 mg/kg”.

## 9 Precision

This Technical Specification has not been tested and validated in an interlaboratory study. Precision values as well as detection limits are therefore not available.

Based on experience, the lowest limit of the nisin A content quantification is 1 mg/kg.

## 10 Test report

The test report shall contain at least the following information:

- a) all information necessary for the complete identification of the sample;
- b) the sampling method used, if known;
- c) the test method used, with reference to this Technical Specification (ISO/TS 27106|IDF/RM 217:2009);
- d) all operational details not specified in this Technical Specification, or regarded as optional, together with details of any incidents which may have influenced the test result(s);
- e) the test result(s) obtained, and, if the repeatability has been checked, the final quoted results obtained.

## Bibliography

- [1] ISO 707|IDF 50, *Milk and milk products — Guidance on sampling*
- [2] ISO 1042, *Laboratory glassware — One-mark volumetric flasks*
- [3] KUZE, N., TAKASHI, M. Analysis of nisin in processed cheese using mass spectrometry coupled with liquid chromatography. *Foods Food Ingredients J. Jpn.* 2003, **208**, 12, pp. 1032-1036
- [4] BERGER, T., BÜTIKOFER, U., MURALT, L., RIEDER, K., RHYN, P., ZÜST, C. Bestimmung von Nisin A in Käse mit LC/MS und LC/MS/MS [Determination of nisin A in cheese with LC-MS and LC-MS-MS]. *Mitt. Lebensm.Hyg.* 2005, **96**, pp. 336-350. Available (2009-10-16) at: [http://www.agroscope.admin.ch/data/publikationen/pub\\_BergerT\\_2005\\_15903.pdf](http://www.agroscope.admin.ch/data/publikationen/pub_BergerT_2005_15903.pdf)

